

360°



VACUETTE[®] QUICKSHIELD Safety Tube Holder

To order, call VWR at 800.932.5000 or visit vwr.com

VACUETTE® QUICKSHIELD Safety Tube Holder

More flexibility with the new **VACUETTE®** QUICKSHIELD Safety Tube Holder

The safety shield of the newly designed holder can now be positioned in any direction. Right-handed and left-handed users can now rotate the safety shield to the position that suits them best.



The rotating safety shield ensures:

- clear visibility of the entire puncture site
- greater flexibility to adapt as effectively as possible to any situation

For daily blood collection procedures

The **VACUETTE®** QUICKSHIELD Safety Tube Holder from Greiner Bio-One is specifically suited for routine blood collection. There is no change to the usual collection technique and one-handed activation of the safety shield is easily accomplished with the aid of a solid surface or with the thumb. This safety device provides the user with ease of handling while protecting the user from exposure to infectious disease.




The ergonomic design of the holder gives fingers a better grip for venipuncture.

The safety shield is pre-assembled directly on the holder. The shield is manually activated upon removal of the needle from the patient's vein.

The **VACUETTE®** QUICKSHIELD Safety Tube Holder and QUICKSHIELD Complete PLUS help prevent needlestick injuries without the need to change collection technique.



VACUETTE® QUICKSHIELD can be used with either **VACUETTE**® Multiple Use Drawing Needles or **VACUETTE**® VISIO PLUS Needles. Both cannulas have a special needle tip bevel that ensures ease of skin puncture. The dot on the protective needle cap is aligned with the needle bevel. This serves as an indication of the bevel position so that the user can adjust the shield to a preferred position without having to remove the protective cap.

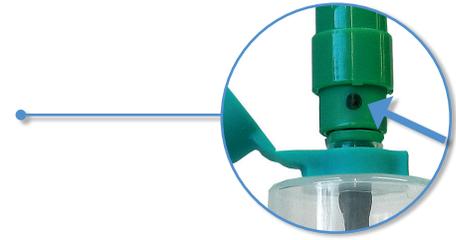

Available as a complete product with pre-assembled **VACUETTE**® VISIO PLUS Needle: **VACUETTE**® QUICKSHIELD Complete PLUS.

Safety Shield
Is pre-assembled directly on the holder and can be positioned in any direction.

VACUETTE® Multiple Use Drawing Needle
Item No. 95057-195

VISIO PLUS
Transparent view window for optical venipuncture control.

VACUETTE® QUICKSHIELD Safety Tube Holder
Item No. 95057-343

VACUETTE® VISIO PLUS Needle Transparent View Window
Item No. 89466-314

VACUETTE® QUICKSHIELD Complete PLUS
The holder comes with pre-attached VISIO PLUS Needle, sterile and packaged in a single blister pack.
Item No. 10789-412, 10789-414, 95057-349, 95057-351

Advantages of the **VACUETTE**® QUICKSHIELD Products:

- Greater flexibility to meet patient needs with a rotating safety shield
- Activated with the aid of a solid surface or thumb
- Convenience through one-handed activation
- Dot on the protective needle cap serves as an orientation aid for position of the needle bevel
- No changes to collection technique
- An audible click indicates to the user that the safety shield has been fully activated
- Can be used with equal ease by both left-handed and right-handed individuals


Instructions for Use *



1 Thread the needle



2 If desired, re-position safety shield

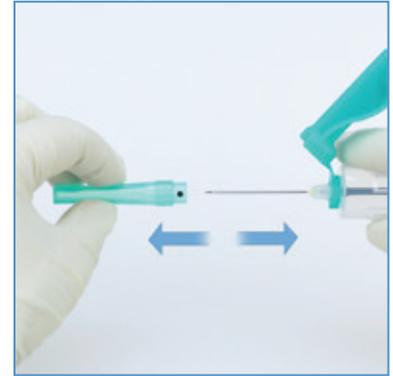

3 Remove protective needle cap



4 Venipuncture



5 Insert the tube



6 Remove the tube



7 Activate the safety mechanism
(with the aid of solid surface or thumb)



8 Audible click indicates activation



9 Disposal

The **VACUETTE**[®] QUICKSHIELD Safety Tube Holder when paired with a **VACUETTE**[®] Needle and the **VACUETTE**[®] QUICKSHIELD Complete PLUS, together with **VACUETTE**[®] Blood Collection Tubes, are suitable as safe, single-use systems for everyday venous blood collection.

* For the complete handling recommendation please refer to the Instructions for Use for each product.

VACUETTE® QUICKSHIELD Safety Tube Holders & Accessories:

					Packaging	
					Inner	Outer
VACUETTE® QUICKSHIELD Safety Tube Holder						
95057-343	450230	VACUETTE® QUICKSHIELD Safety Tube Holder			50	1,000
VACUETTE® Multiple Use Drawing Needles (sterile)						
95057-195	450071	VACUETTE® Multiple Use Drawing Needle, black	0.7 x 25mm	22G x 1"	100	2,000
95057-197	450072	VACUETTE® Multiple Use Drawing Needle, green	0.8 x 25mm	21G x 1"	100	2,000
95057-199	450073	VACUETTE® Multiple Use Drawing Needle, yellow	0.9 x 25mm	20G x 1"	100	2,000
95057-201	450075	VACUETTE® Multiple Use Drawing Needle, black	0.7 x 38mm	22G x 1 1/2"	100	2,000
95057-203	450076	VACUETTE® Multiple Use Drawing Needle, green	0.8 x 38mm	21G x 1 1/2"	100	2,000
95057-205	450077	VACUETTE® Multiple Use Drawing Needle, yellow	0.9 x 38mm	20G x 1 1/2"	100	2,000
VACUETTE® VISIO PLUS Needles (sterile)						
89166-314	450040	VACUETTE® VISIO PLUS Needle, green	0.8 x 38mm	21G x 1 1/2"	100	2,000
89166-316	450041	VACUETTE® VISIO PLUS Needle, black	0.7 x 38mm	22G x 1 1/2"	100	2,000
89166-318	450042	VACUETTE® VISIO PLUS Needle, green	0.8 x 25mm	21G x 1"	100	2,000
89166-320	450043	VACUETTE® VISIO PLUS Needle, black	0.7 x 25mm	22G x 1"	100	2,000

These devices have no components made of natural dry rubber.

Your **Power** for Health


greiner bio-one

360°



980959 rev.01, 03.2015 e - VWR

VWR  **1.800.932.5000**
We Enable Science **vwr.com**



Prices and product details are current when published; subject to change without notice. | Certain products may be limited by federal, state, provincial, or local regulations. | VWR makes no claims or warranties concerning sustainable/green products. Any claims concerning sustainable/green products are the sole claims of the manufacturer and not those of VWR International, LLC. All prices are in US dollars unless otherwise noted. Offers valid in US, void where prohibited by law or company policy, while supplies last. | VWR, the VWR logo and variations on the foregoing are registered (®) or unregistered trademarks and service marks, of VWR International, LLC and its related companies. All other marks referenced are registered by their respective owner(s). | Visit vwr.com to view our privacy policy, trademark owners and additional disclaimers. ©2015 VWR International, LLC. All rights reserved.