

Certificate of Analysis

BamHI Restriction Enzyme:

Part No.	Size (units)
R602A	2,500
R602B	12,500
R402A	(High Conc.) 12,500
R402B	(High Conc.) 50,000

Buffer E (R005A): The Restriction Enzyme 10X Buffer E supplied with this enzyme has a composition of 60mM Tris-HCl (pH 7.5), 1M NaCl, 60mM MgCl₂ and 10mM DTT at 37°C.

BSA in Restriction Enzyme Reactions: All Promega restriction enzyme activity assays are performed with the addition of Acetylated BSA to a final concentration of 0.1mg/ml. BSA has been shown to enhance the activity of many restriction enzymes. Promega supplies a tube of Acetylated BSA (Part# R396D, E or F) with all restriction enzymes and recommends the addition of BSA to a final concentration of 0.1mg/ml to all restriction digestions.

Assay Information

Overdigests (OD): The overdigestion assay is performed to demonstrate the absence of nonspecific endonucleases and exonucleases. Increasing amounts of enzyme are added to a series of tubes containing substrate DNA. After a 16-hour incubation under appropriate conditions, the maximum number of units giving a clear, sharp, normal banding pattern is determined by agarose gel/ethidium bromide electrophoresis. The minimum number of units of BamHI required to meet these criteria using 1µg of lambda DNA as the substrate and the recommended digestion conditions (50µl reaction, Buffer E, 37°C) is 20 units.

Cut:Ligation:Recut (C:L:R): The ligation assay is used to determine the functional purity of the DNA after restriction enzyme digestion. Substrate DNA is completely digested with a fourfold excess of the restriction enzyme in the appropriate assay buffer, ligated with T4 DNA Ligase and recut with the same restriction enzyme. Cut, ligated and recut DNAs are analyzed by agarose gel/ethidium bromide electrophoresis.

Product Information

Buffer Activity:

Percent (%) activity in 4-CORE® Buffer System

A	B	C	D	E	MC
75-100	75-100	75-100	50-75	100	75-100

Key: A-D 4-CORE® Buffer System
MC MULTI-CORE™ Buffer

Heat Inactivation:

	Heat Inactivated
BamHI	+
Key:	+ 95% inactivation by incubation at 65°C for 15 minutes
	- no inactivation

MULTI-CORE™ Buffer (R999A): The MULTI-CORE™ 10X Buffer, which is supplied with enzymes that exhibit >25% activity in this buffer, has a composition of 250mM Tris acetate (pH 7.8), 1M potassium acetate, 100mM magnesium acetate and 10mM DTT. BamHI is 75-100% active in MULTI-CORE™ Buffer containing 0.1mg/ml Acetylated BSA.

Recognition Sequence : 5'...GGATC C...3
3'...C CTAG G...5

Enzyme Storage Buffer: BamHI is supplied in 10mM Tris-HCl (pH 7.4), 300mM KCl, 5mM MgCl₂, 0.1mM EDTA, 1mM DTT, 0.5mg/ml BSA and 50% glycerol.

Source: *Bacillus amyloliquefaciens* H. (1).

Unit Definition: One unit is defined as the amount of enzyme required to completely digest 1µg of lambda DNA in one hour at 37°C in 50µl assay buffer containing Acetylated BSA added to a final concentration of 0.1mg/ml.

Optimal Assay Temperature: 37°C.

Storage Temperature: See Product Information Label for storage recommendations and expiration date. Avoid frequent temperature changes.

A normal banding pattern indicates intact 5' and 3' termini as well as the absence of contaminating nucleases. The minimum C:L:R ratio for this enzyme is 100:90:90.

Blue/White Cloning Qualified: This assay is performed to determine the incidence of false positives (2). A pGEM®-Z series Vector containing a unique BamHI site is digested with 25 units of BamHI per 5µg of DNA for 2 hours at 37°C. Linearized plasmid is purified using Promega's Wizard® Plus Minipreps DNA Purification System and ligated at 25°C for 30 minutes. Competent JM109 cells are transformed with ligated plasmids and plated on X-Gal/IPTG/Amp₁₀₀ plates. A minimum of 400 colonies are counted. White colonies result from transformation with ligated plasmids with damaged ends. These white colonies represent the number of false positives expected in a typical cloning experiment. Enzymes that generate overhangs must produce fewer than 2% white colonies, and blunt-cutting enzymes must produce fewer than 5% white colonies.

+/- partial inactivation

Frequency of Cutting:

	1	Ad-2	ϕX174	pUC18	M13mp18	pBR322
5	3	0	1	1	1	1

Methylation:


	dam	dcm	CpG	CpNpG
BamHI	i	i	i	s(ol)

Key: s = sensitive to this methylation
i = insensitive to this methylation
s(ol) = overlapping (sensitive when restriction site overlaps methylation sequence)
n/a = information not available

References

- Roberts, R.J. *et al.* (1977) *Nature* 265, 82.
- Hung, L. *et al.* (1991) *Promega Notes* 33, 12.

Signed by:

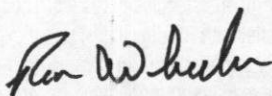

R. Wheeler, Quality Assurance

BamHI

REF R6021

LOT 0000305939

-30°C -10°C

2020-08-30

Conc: 10U/µl
2,500U

Dispensed Lot#: 0000290154

For In Vitro Research Use Only

Country of Origin: USA

Promega Corporation
2800 Woods Hollow Road
Madison, WI 53711-5399 USA



ADR6021 00003059394

PEEL
HERE



Promega Corporation

2800 Woods Hollow Road	
Madison, WI 53711-5399 USA	
Telephone	608-274-4330
Toll Free	800-356-9526
Fax	608-277-2516
Internet	www.promega.com

PRODUCT USE LIMITATIONS, WARRANTY, DISCLAIMER

Promega manufactures products for a number of intended uses. Please refer to the product label for the intended use statements for specific products. Promega products contain chemicals which may be harmful if misused. Due care should be exercised with all Promega products to prevent direct human contact.

Each Promega product is shipped with documentation stating specifications and other technical information. Promega products are warranted to meet or exceed the stated specifications. Promega's sole obligation and the customer's sole remedy is limited to replacement of products free of charge in the event products fail to perform as warranted. Promega makes no other warranty of any kind whatsoever, and SPECIFICALLY DISCLAIMS AND EXCLUDES ALL OTHER WARRANTIES OF ANY KIND OR NATURE WHATSOEVER, DIRECTLY OR INDIRECTLY, EXPRESS OR IMPLIED, INCLUDING, WITHOUT LIMITATION, AS TO THE SUITABILITY, PRODUCTIVITY, DURABILITY, FITNESS FOR A PARTICULAR PURPOSE OR USE, MERCHANTABILITY, CONDITION, OR ANY OTHER MATTER WITH RESPECT TO PROMEGA PRODUCTS. In no event shall Promega be liable for claims for any other damages, whether direct, incidental, foreseeable, consequential, or special (including but not limited to loss of use, revenue or profit), whether based upon warranty, contract, tort (including negligence) or strict liability arising in connection with the sale or the failure of Promega products to perform in accordance with the stated specifications.

© 1996, 1998-2000, 2004, 2005, 2008, 2013, 2016
Promega Corporation. All Rights Reserved.

MULTI-CORE is a trademark of Promega Corporation. 4-CORE, pGEM and Wizard are registered trademarks of Promega Corporation.

Products may be covered by pending or issued patents or may have certain limitations. Please visit our Web site for more information.

All specifications are subject to change without prior notice.

Product claims are subject to change. Please contact Promega Technical Services or access the Promega online catalog for the most up-to-date information on Promega products.

Part# 9PIR602
Printed in USA. Revised 8/16.

Introduction

Restriction enzymes, also referred to as restriction endonucleases, are enzymes which recognize short, specific (often palindromic) DNA sequences. They cleave double-stranded DNA (dsDNA) at specific sites within or adjacent to their recognition sequences. Most restriction enzymes (REs) will not cut DNA that is **methylated** on one or both strands of their **recognition site**, although some require substrate methylation.

Each restriction enzyme has specific requirements to achieve optimal activity. Ideal storage and assay conditions favor the most activity and highest fidelity in a particular enzyme's function. Conditions such as temperature, pH, enzyme cofactor(s), salt composition and ionic strength affect enzyme activity and stability. Two buffers usually accompany each of Promega's restriction enzymes. One buffer is the optimal reaction buffer which may be from the **4-CORE® System** (Reaction Buffers A, B, C, D) or one of the other optimal buffers (Reaction Buffers E–L), and the other is the **MULTI-CORE™ Buffer**. The supplied optimal buffer always yields 100% activity for the enzyme it accompanies, and serves as the specific reaction buffer for individual digests with that enzyme.

The **MULTI-CORE™ Buffer**, which is designed for broad compatibility with many REs, is provided with enzymes that have 25% or greater activity in this buffer. The **MULTI-CORE™ Buffer** is useful for multiple digests because it generally yields more activity for more enzyme combinations than any of the other buffers, but sometimes with a compromise in activity. Multiple digests using REs with significantly different buffer requirements may require a sequential reaction with the addition of RE buffer or salt before the second enzyme is used.

DNA Substrate Considerations

DNA substrates commonly used for restriction enzyme digestion include DNA from bacteriophage lambda, bacterial plasmid DNA and genomic DNA. Lambda DNA is a linear DNA form that is an industry standard for measuring and expressing unit activity for many restriction enzymes. Compared to linear DNA, intact supercoiled plasmid DNA (and DNAs with a large number of the target restriction site) require more units of enzyme (two- to tenfold) per microgram than the DNA used in the enzyme's activity assay.

DNA Substrate Considerations

PCR products and oligonucleotides are relatively small compared with DNA used for defining RE units. Therefore, when using these substrates in a restriction digest, it is essential to take into consideration the molar concentration of enzyme recognition sites and not just the mass of DNA. Also, some REs require flanking bases surrounding the core RE recognition site. This is problematic when it is necessary to cut an oligonucleotide or a fragment of DNA with an RE site near its end. When PCR cloning strategies include the use of primers containing an RE site, care is necessary in designing the primer with adequate DNA surrounding the core RE recognition sequence.

In addition to the form and original source of the DNA, the purity is another factor that must be considered. Depending on the purification method and the handling of the DNA, it may contain varying amounts of contaminants that affect restriction enzyme digestion and analysis. Contaminants may include other types of DNA, nucleases, salts and inhibitors of restriction enzymes. The effect of a contaminant on an RE digest is generally dose-dependent: i.e., the inhibitory effects will increase with the volume of DNA added to the restriction enzyme reaction. Relatively pure DNA is required for efficient restriction enzyme digestion. Contaminating nucleases are usually activated only after the addition of salts (e.g., restriction enzyme buffer) to the DNA solution. Therefore, appropriate control reactions should always be run in parallel with the restriction digest. Buffer solutions containing EDTA in low concentrations (1mM) are often used to protect DNA from nuclease degradation during storage,

but the EDTA can interfere with restriction enzyme digestion if the final concentration of EDTA in the reaction is too high. This situation usually results when the concentration of the substrate DNA is low and it is necessary to use a large volume of DNA in the digest. In such cases, it is best to concentrate the DNA (e.g., by ethanol precipitation). The organic solvents, salts, detergents and chelating agents that are sometimes used during the purification of DNA can also interfere with restriction enzyme activity if they carry over into the final DNA solution. Dialysis and/or ethanol precipitation with 2.5M ammonium acetate (final concentration before adding ethanol) followed by drying and resuspension can remove many of these substances. While relatively pure DNA is required for efficient restriction enzyme digestion, addition of acetylated BSA to a final concentration of 0.1mg/ml can sometimes improve the quality and efficiency of enzyme assays containing impure DNA and we recommend that it be included in all digests.

Enzyme Storage, Handling and Use

Maintain the sterility of reagents used in the RE digest as well as any tools (e.g., tubes, pipette tips) used with those reagents. Restriction enzymes should be stored in a non-frost-free freezer, except for a brief period during use, when they should be kept on ice. The restriction enzyme is usually the last reagent added to a reaction, to ensure that it is not exposed to extreme conditions. When many similar digests are being prepared, it may be convenient to create pre-mixes of common reagents.

Before assembling the restriction digest, thoroughly mix each component to be added to the reaction and then centrifuge the tubes of reagents briefly to collect the contents in the bottom of the tube. The reaction components should also be mixed after addition of the enzyme to the digest. While high salt buffers and glycerol-containing reagents are difficult to mix, all solutions containing restriction enzymes must be mixed gently to avoid inactivating the enzyme.

Enzyme Storage, Handling and Use

Setting up a Restriction Enzyme Digest

An analytical scale restriction enzyme digest is usually performed in a volume of 20µl on 0.2–1.5µg of substrate DNA, using a two- to tenfold excess of enzyme over DNA. If an unusually large volume of DNA or enzyme is used, aberrant results may occur and may or may not be readily recognized. The following is an example of a typical RE digest. In a sterile tube, assemble in order:

sterile, deionized water	16.3µl
RE 10X Buffer	2µl
Acetylated BSA, 10µg/µl	0.2µl
DNA, 1µg/µl	1µl

Mix by pipetting, then add:

Restriction Enzyme, 10u/µl	0.5µl
final volume	20µl

Mix gently by pipetting, close the tube and centrifuge for a few seconds in a microcentrifuge. Incubate at the optimum temperature for 1–4 hours.

Add 4µl of 6X loading buffer and proceed to gel analysis. Note that overnight digests are usually unnecessary and may result in degradation of the DNA.

Setting up a Restriction Enzyme Digest

Experimental Controls

Experimental controls are necessary to identify, understand and explain problems or inconsistencies in results. The following controls are commonly used in parallel with RE digests: (i) uncut experimental DNA, (ii) digest of commercially supplied control DNA, (iii) no-enzyme "mock" digest, (iv) 1 or 2 different size markers in more than one lane per gel

Part# 9PIR602
Printed in USA. Revised 8/16.